

SURIMEX

Surgical Import Export Distributor

Safety Infusion Set

›EZ Huber™

The EZ Huber™ is a safety needle for puncturing implanted port systems.

A safety mechanism reduces the risk of being injured by the needle and protects users from infection.

www.surimex.com.au



The **EZ Huber™** is a safety needle for puncturing implanted port systems. A safety mechanism reduces the risk of being injured by the needle and protects users from infection.

View



(1) prior use

(2) after use

Benefits

Patented safety mechanism

The patented EZ Huber™ safety mechanism reduces the risk of needlestick injuries and provides protection against infection. The safety mechanism is activated automatically when the cannula is extracted.

An audible click indicates that the cannula is safely locked inside the housing. The EZ Huber™ meets the requirements of the TRBA 250 for security products.

(TRBA 250: German technical regulations governing biological working materials in health care and welfare work.)

Unique protective sheath

The unique, patented protective sheath completely covers the cannula when it is withdrawn, providing protection against contact with bodily fluids and germs. The protection sheath effectively prevents cross contaminations.

Biocompatible material

The use of extra biocompatible materials - containing no latex and DEHP- affords the patient additional protection. Infusion of aggressive medication without material damage is still possible.


Details

- ▶ The EZ Huber™ is suitable for high-pressure injection procedures with pressures up to 21 bar / 300 psi , e.g. for injection of contrast media during CT-Scans (only in combination with an appropriate port system, such as T-Port Contrast and Jet Port® Plus II Contrast developed by pfm medical ag). The EZ Huber™ is universally recognised as high-pressure, by its colour purple.
- ▶ The wings of the EZ Huber™ provide a good grip to ensure convenient puncturing procedures and easy removal of the needle.
- ▶ The configuration of the wings and the transparent base ensures a visual inspection of the injection site to easily detect any signs of infection.
- ▶ The conical shape of the foam padding prevents pressure marks or irritation of the patient's skin.
- ▶ The very low profile of EZ Huber provides a safe and reliable fixation and increases the patient's comfort.
- ▶ The adjustable base is able to level differences between the base and the skin and provides an increased stability during use.
- ▶ The EZ Huber™'s simple handling does not require any changes in the usual routine with respect to the use of port needles.
- ▶ The EZ Huber™ with Y-site* allows additional medication to be administered without interrupting an infusion in progress.




* Only available upon request

Handling

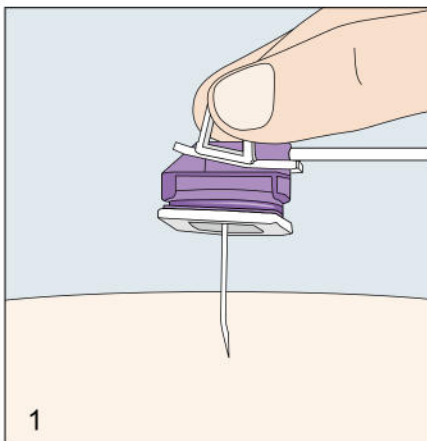


Insert

- ▶ Remove needle cap.
- ▶ Hold EZ Huber™ at the wings and insert needle into the skin over the port.
- ▶ Ensure needle has reached bottom of port.

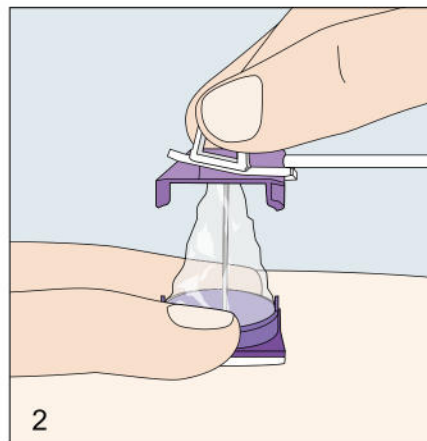

Remove

- ▶ Stabilise the port by placing two fingers on the base of the EZ Huber™.
- ▶ Withdraw the needle by pulling up on the wings of the EZ Huber™.

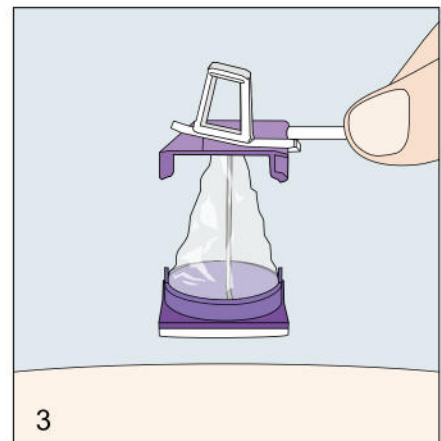

Dispose

- ▶ Dispose EZ Huber™ as per institutional protocol.

Ordering Information

Ez Safety Huber Needle (Extension 180mm)

REF	DESCRIPTION	Size	Colour	PU	Stock Availability*
SHG19-75	Ez Safety Huber Needle	19G x 19mm	■ Brown	25	Available upon request
SHG19-100	Ez Safety Huber Needle	19G x 25mm	■ Brown	25	Available upon request
SHG19-150	Ez Safety Huber Needle	19G x 38mm	■ Brown	25	Available upon request
SHG20-75	Ez Safety Huber Needle	20G x 19mm	■ Yellow	25	Stocked Item
SHG20-100	Ez Safety Huber Needle	20G x 25mm	■ Yellow	25	Stocked Item
SHG20-150	Ez Safety Huber Needle	20G x 38mm	■ Yellow	25	Stocked Item
SHG20-75Y	Ez Safety Huber Needle with Y-Site Connector	20G x 19mm	■ Yellow	25	Stocked Item
SHG20-100Y	Ez Safety Huber Needle with Y-Site Connector	20G x 25mm	■ Yellow	25	Stocked Item
SHG22-75	Ez Safety Huber Needle	22G x 19mm	■ Black	25	Available upon request
SHG22-100	Ez Safety Huber Needle	22G x 25mm	■ Black	25	Available upon request
SHG22-150	Ez Safety Huber Needle	22G x 38mm	■ Black	25	Available upon request

*For current stock availability, please check our website: www.surimex.com.au or contact our Head Office on (03) 8833 3804

Contact person

Do you have any questions? Our Customer Sales & Service Support team will be glad to assist you.

Products Distributed in Australia and New Zealand by
www.surimex.com.au



15 Redland Drive
 Mitcham VIC 3132
 Australia
 T (613) 8833 3804
 E sales@surimex.com.au