


SURIMEX

Drains

PVC, Silicone and Silicone Fluted

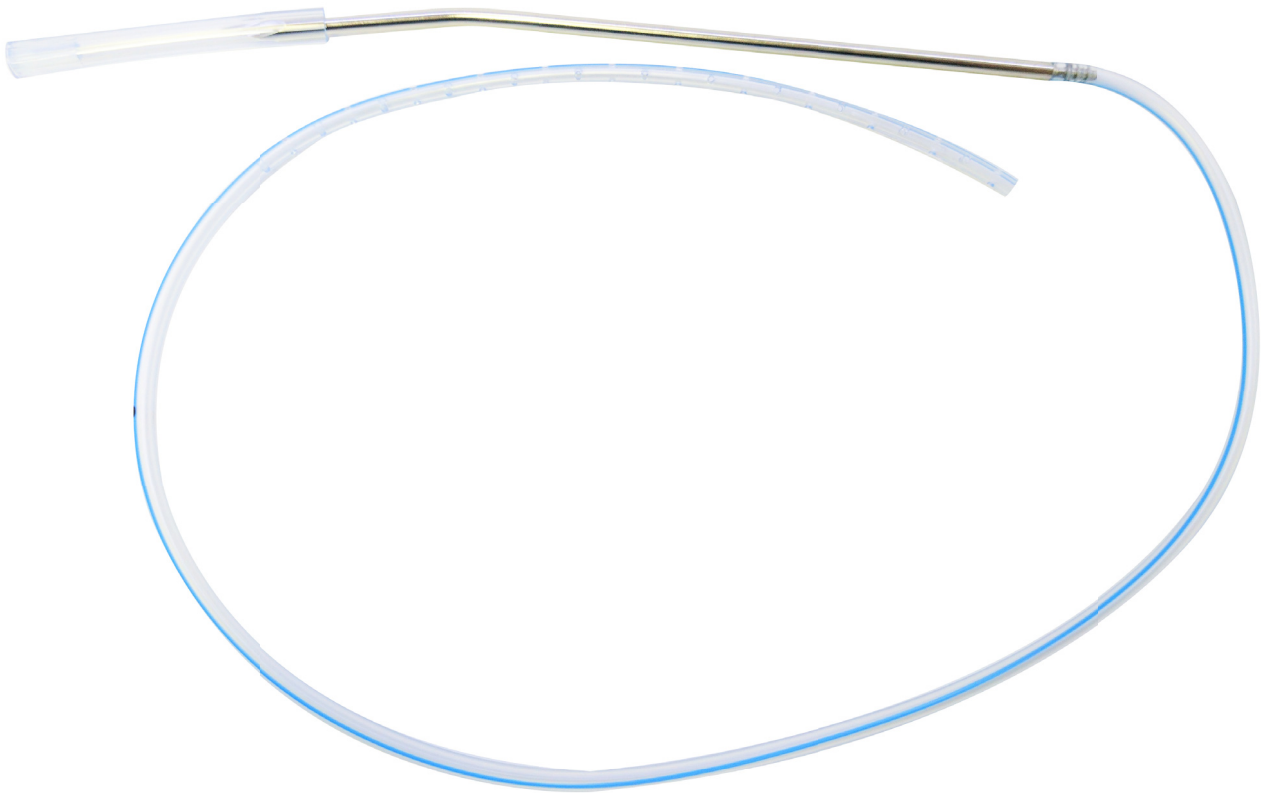


Drains PVC, Silicone and Silicone Fluted

Our drain range is for post operative use with high, low or gravity drainage systems

PVC drains: are cross perforated, suitable for use with high / low vacuum wound drainage systems.


Silicone drains: are cross perforated, suitable for use with high / low vacuum wound drainage systems.

Silicone fluted drains: are a 4 channelled drain, (generally used for sensitive areas), suitable for use with high / low vacuum or gravity wound drainage systems.

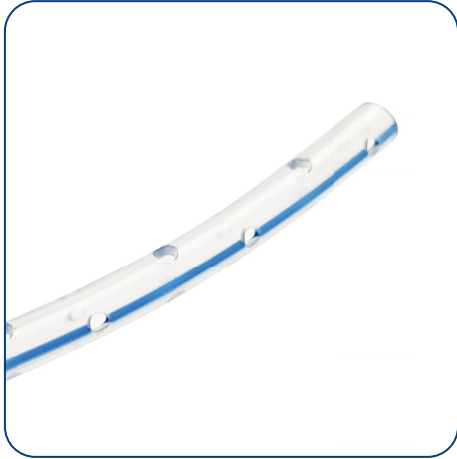


Cross perforation

- Perforations along distal end

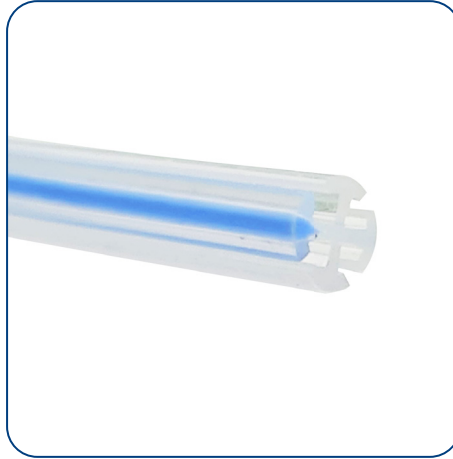

Silicone fluted drain channels

- Kink resistant due to fluted channels
- Reduces tissue ingrowth and blockage

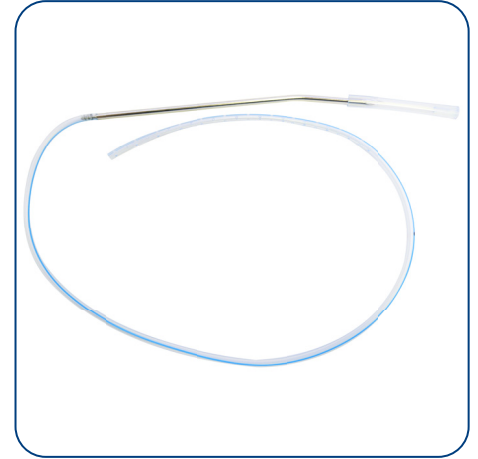

Pre-connected trocar

- Stainless steel, sharp, bur free
- Angled for ease of use



Radiopaque lines

- X-ray visualisation



Silicone drain connector

- Secure, tapered design
- Allowing for optimal flow

Item no.	Description	Drain FG	Perforation length	Drain length	radiopaque lines	with trocar	with connector	PU	Stock availability
220021-P	PVC Drain (drain only)	10	15cm	50cm	•			150	stocked item
221510	PVC Drain Set	10	15cm	50cm	•	•		10	stocked item
221514	PVC Drain Set	14	20cm	50cm	•	•		10	stocked item
221610	Silicone Drain Set	10	15cm	50cm	•	•	•	10	stocked item
S33501-P	Silicone Fluted Drain Set	10	20cm	110cm	•	•	•	10	stocked item
S33502-P	Silicone Fluted Drain Set	15	20cm	120cm	•	•	•	10	stocked item

Do you have any questions? Our customer service & sales support will be glad to assist you.

Products distributed in Australia and New Zealand by



15 Redland Drive
Mitcham VIC 3132
Australia

T (613) 8833 3804
E sales@surimex.com.au
www.surimex.com.au